

Product Codes

Sizing Information

Size	Length of sling (mm)	Code
Small	340 - 380	0663
Medium	370 - 430	0665
Large	420 - 480	0667
XLarge	480+	0668



BEA SHOULDER SLING FITTING INSTRUCTIONS

Any product supplied by Beagle Orthopaedic is intended for single patient use only. Re-use of any product will invalidate the CE mark and any warranty.




Beagle Orthopaedic


Prospect Building, Duttons Way, Shadsworth Business Park, Blackburn, BB1 2QR
E: info@beagleortho.com W: beagleorthopaedic.com T: 01254 268 788

INDICATIONS FOR USE


» Immobilisation following injury or surgery

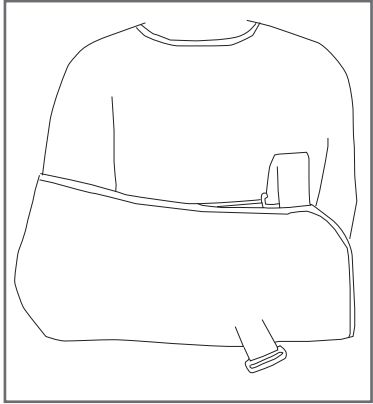


1 Place the injured arm into the sling, ensuring the elbow and forearm are supported securely to avoid discomfort

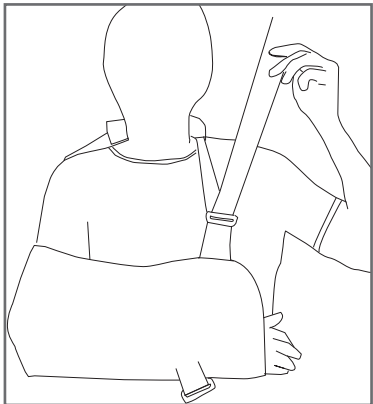


2 Run the shoulder strap up from the injured elbow, round the neck and over the uninjured shoulder. Feed it through the D-ring on the outside of the sling and secure with the hook & loop fastening

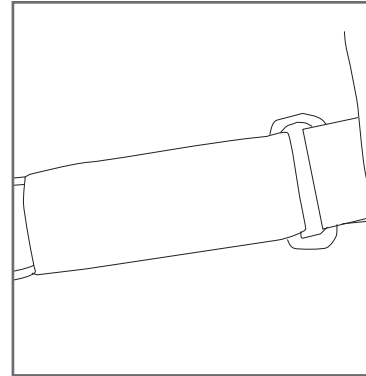

3 Attach the immobilisation strap on to the side of the sling

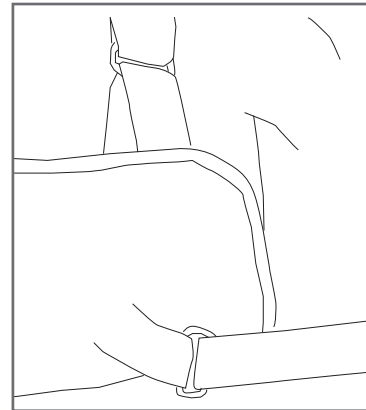

4 Feed the immobilisation strap through the D-ring (as shown) and secure with the hook & loop fastening.

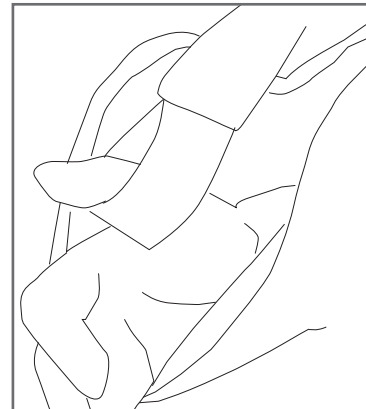

5 Soft thumb loop maintains optimal position and supports hand